

# AH Plus.<sup>®</sup> The Superior Sealer.



## *AH Plus Root Canal Sealer Is A Resin-Based Formula With:*

- Excellent Radiopacity
- Low Shrinkage
- Low Solubility
- Outstanding Flow
- A 1:1 Paste-To-Paste System For Fast, Easy Mixing And Less Waste

### *AH Plus Is Safe And Biologically Compatible.*

AH Plus's patented formula uses new types of amines that result in excellent biocompatibility. It is chemically inert after curing, making it a viable solution to postoperative problems such as sensitivity and periapical inflammation.

### *Designed For Excellent Flow And Low Shrinkage And Solubility.*

AH Plus has outstanding physical properties. Its excellent flow characteristics mean AH Plus can reach and seal lateral canals. Low solubility makes it resistant to break down over time. And to reduce gapping between the sealer and the canal wall, AH Plus exhibits minimal shrinkage.

### *Whatever Your Technique, The Seal Remains The Same.*

AH Plus has been formulated to work well with any canal sealing technique—from the traditional to the most innovative. AH Plus works especially well with warm gutta-percha methods because it has excellent flow characteristics when heated.


### *Getting The Right Mix Is As Easy As 1:1.*

With most sealers, there is an art to getting the right mix to achieve the proper consistency. With AH Plus, there's a science. Its paste/paste system with 1:1 mixing ratio gives you more precise measuring prior to mixing. This eliminates guesswork and provides better dosage control and less waste.

**DENTSPLY  
MAILLEFER**

*In Endodontics, Quality Has A Name.*



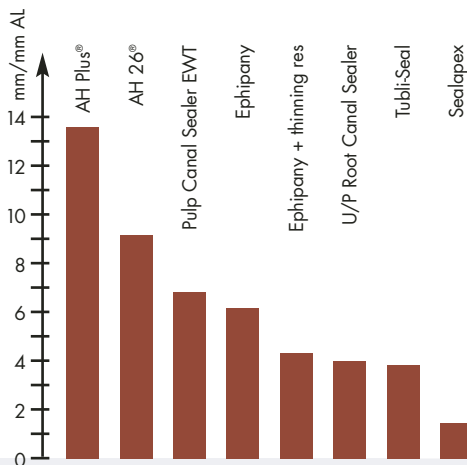

AH Plus performed best in  
**Tests Measuring Radiopacity.**

## AH Plus Gives You Optimum Radiopacity.

AH Plus contains dense radiopaque fillers to ensure exceptionally good radiopacity even when applied in thin layers. In fact, in laboratory testing, AH Plus showed the highest radiopacity when tested against the five other leading brands.



Radiopacity	13.6 mm/mm Al
Working time	4 hours
Setting time	8 hours
Flow	36 mm
Film thickness	26 µm
Shrinkage	1.76%
Solubility	0.31%

**Technical data**

## Easy To Work With From Start To Finish.

With a working time of four hours and a setting time of eight hours, AH Plus is easy to use. And its thermoplastic properties make it easier to remove, whether you're using chemicals or a drill.

## Achieve An Improved Seal With The Superior Sealer.

For more information about AH Plus or for the nearest Maillefer distributor, contact us.

Features:	Benefits:
Paste/paste system	Fast, easy mixing and dosage control
Excellent radiopacity	Easy detection
Outstanding dimensional stability and low solubility	Durable, safe seal
Outstanding flow characteristics	Flows into and seals lateral canals
Biologically inert	Biocompatible
New amines in formula (silver free)	No discoloring of teeth
Thermoplastic ability	Easier to remove
Four-hour working time	No waste of material on pad
Eight-hour setting time	Allows for corrections of fillings
Compatible with any filling technique	Versatility



800-924-7393 [www.maillefer.dentsply.com](http://www.maillefer.dentsply.com) 800-924-7389 (fax)

5100 East Skelly Drive, Suite 300, Tulsa, OK 74135