

All the Characteristics You Know and Trust.

While AH Plus Jet sealer delivers added efficiency and precision, its preloaded sealer remains true to its proven, patented formula. You get all the properties endodontists have come to expect – starting with excellent biocompatibility, making it a viable solution to postoperative problems such as sensitivity and periapical inflammation. In addition, its low solubility makes it resistant to breaking down over time. And its minimal shrinkage reduces gapping between the sealer and the canal wall. Should the need for removal arise, the sealer's thermoplastic properties make the task easier, whether using chemicals or a drill.

AH Plus Jet Root Canal Sealer

Description	Item#
AH Plus® Jet™ Starter Kit	667005
AH Plus® Jet™ Refill Tips	667006


Technical Data

Radiopacity	13.6 mm/mm Al
Working time	4 hours
Setting time	8 hours
Flow	36 mm
Film thickness	26 µm
Shrinkage	1.76%
Solubility	0.31%

The Smart Sealer. AH Plus Jet Root Canal Sealer.

For more information, call your authorized Maillefer distributor.

800-924-7393 www.maillefer.dentsply.com 800-924-7389 (fax)
5100 East Skelly Drive, Suite 300, Tulsa, OK 74135

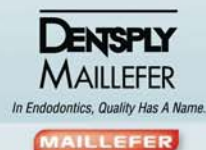

AH Plus® Jet™ Root Canal Sealer



Self-Mixing Efficiency.

The Smart Sealer: AH Plus Jet Root Canal Sealer.

- Ready-to-use efficiency with no manual mixing
- Controlled, homogeneous, automatic mixing of both pastes
- Better working ergonomics for direct, precise placement
- Simple application, no assistance needed
- Intra Oral Tip adjustable to individual anatomic conditions
- Clean application
- Excellent radiopacity
- Low shrinkage
- Low solubility
- Outstanding flow characteristics



AH Plus[®] Jet[™]

Root Canal Sealer

Same Great Seal. Even Greater Efficiency.


Achieving an excellent seal with the time-tested formula endodontists trust is more efficient with AH Plus Jet sealer. Its innovative delivery system eliminates the need for manual mixing prior to use, while enabling direct and precise placement into the canal or on a traditional mixing pad.

AH Plus Jet Sealer Does The Mixing For You.

Utilizing the proven, resin-based AH Plus sealer formula, AH Plus Jet sealer features an innovative double-barrel syringe that automatically and precisely mixes the paste-to-paste formula in the necessary 1:1 ratio. No manual mixing or assistance is required. Just depress the syringe's plunger for a controlled, homogenous mixture – and the confidence of knowing it's exactly right.

Easy to Use, Whatever Your Technique.

With a working time of four hours and a setting time of eight hours, AH Plus Jet sealer is designed to be easy to work with and effective. Its excellent flow characteristics allow it to reach and seal lateral canals, making it especially suited to warm gutta-percha methods.

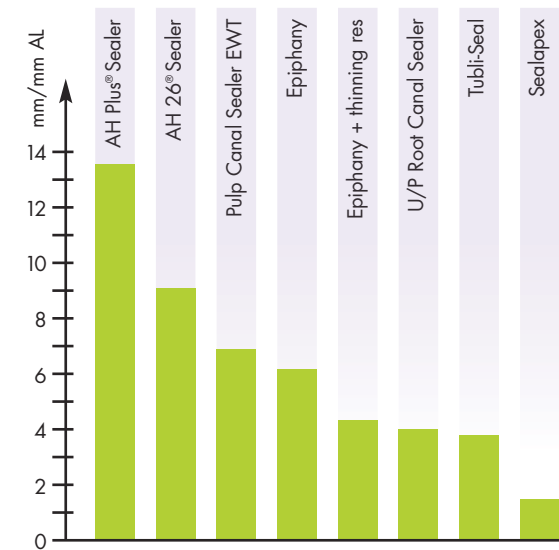

AH Plus Jet sealer performed best in Tests Measuring Radiopacity.*

Optimum Radiopacity.

AH Plus Jet sealer contains dense radiopaque fillers to ensure exceptionally good radiopacity. In fact, in laboratory testing, the sealer showed the highest radiopacity when tested against other leading brands.



*Study performed by DENTSPLY DeTrey, Konstanz, Germany.